

Phospho-p44/42 MAPK (Erk1/2) (Thr202/Tyr204) (A11) rabbit mAb FITC conjugate

www.abwizbio.com
Support: info@abwizbio.com
Order: sales@abwizbio.com


Catalog: #1113

Store at: 2-8°C

For Research Use Only. Not For Use In Diagnostic Procedures.

Applications	Detection	Clonality	Isotype
Flow Cytometry	N/A	Monoclonal	Rabbit IgGk

Format:	FITC
Cross Reactivity:	Predicted to work with mouse, rat, and other homologues.
Formulation:	PBS, 0.09% NaN ₃ , 0.2% BSA
Preparation:	Protein A
Reactivity:	Human, Mouse
Recommended Usage:	For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells or 5 µL per 100 µL of staining volume. It is recommended that the reagent be titrated for optimal performance for each application.
Immunogen:	A synthetic phospho-peptide corresponding to residues surrounding Thr202/Tyr204 of human phospho Erk1/2.
Description:	Human Erk1 and Erk2 Ser/Thr kinases share 84% sequence identity and nearly all functions. These MAP kinases are activated in response to mitogens and growth factors as part of the Ras-Raf-MEK-ERK signal transduction cascade(1-3). This pathway regulates cell survival, differentiation, adhesion, cell cycle progression, and many other cellular processes. Upon phosphorylation, Erk1/2 translocate to the nucleus to activate transcription factors including c-Fos, Elk1, Ets1, and SP-1 (4,5). There are more than 175 known cytoplasmic and nuclear substrates of Erk1/2. The Erk1/2 cascade is upregulated in many human cancers, even when oncogenic mutations are not found. Multiple small-molecule inhibitors of Erk1/2 have been developed, including ones targeting the ATP-binding site either competitively or irreversibly (6).
References:	<ol style="list-style-type: none">1. Blagoev B, et al., 2003, Nat Biotechnol, 21:315-318.2. Thelemann A. et al, 2005, Mol Cell Proteomics 4:356-376.3. Morandell S., et al., 2008, Proteomics 8:4383-4401.4. Ramos J.W., 2008, Biochem Cell Biol. 40:2707-2719.5. Nakano H., et al., 1998, Proc Natl Acad Sci U S A. 104:19837-19842.6. Roskoski Jr R. 2012, Pharmacol Res. 66: 105-143.

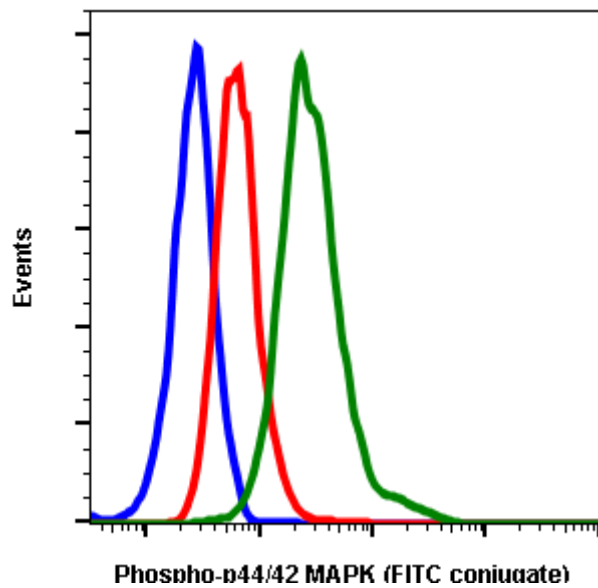

Flow cytometric analysis of Jurkat cells secondary antibody only negative control (blue) or treated with U0126 (red) or treated with TPA (green) using Phospho-ERK1/2 (Thr202/Tyr204) A11 FITC antibody ERK12T202Y204-A11. Cat. #1113.