

Phospho-Akt1 (Ser473) (C7) rabbit mAb

www.abwizbio.com

Support: info@abwizbio.com

Order: sales@abwizbio.com

Catalog: #1131

Store at: -20°C


For Research Use Only. Not For Use In Diagnostic Procedures.

Applications	Detection	Clonality	Isotype
Flow Cytometry	Anti-Rabbit IgG	Monoclonal	Rabbit IgGk


Format:	Unconjugated
Cross Reactivity:	Predicted to work with mouse, rat and other homologues.
Formulation:	1X PBS, 0.02% NaN ₃ , 50% Glycerol, 0.1% BSA
Preparation:	Protein A+G
Reactivity:	Human, Mouse
Recommended Usage:	1µg/mL - 0.001µg/mL. It is recommended that the reagent be titrated for optimal performance for each application. See product image legends for additional information.
Immunogen:	A synthetic phospho-peptide corresponding to residues surrounding Ser473 of human phospho Akt1
Description:	<p>Akt also known as PKB (Protein kinase B) or RAC-PK (Related to the A and C kinases) is a serine/threonine kinases that contains a pleckstrin homology (PH) domain. This protein kinase is activated by insulin and various growth and survival factors to function in a wortmannin sensitive pathway involving PI3 kinase. Akt is activated by phospholipid binding and activation loop phosphorylation at Thr308 by PDK1 and by phosphorylation within the carboxy terminus at Ser473. Phospho-Akt promotes cell survival by inhibiting apoptosis. Specifically, phospho-Akt1 has been shown to phosphorylate Bad, a member of the Bcl-2 family that promotes cell death. This phosphorylation results in the inactivation of the proapoptotic function of Bad. The Akt/phospho Akt molecule is thus considered to link extracellular survival signals (growth factors) with the apoptotic machinery (BAD). Akt is also a key mediator of the metabolic effects of insulin. Additionally, Akt has been referred to as an oncogene because it has increased activity in a number of tumors. This antibody recognizes phospho Akt phosphorylated at Ser473. This phosphorylation site is shared by all three isoforms of phospho Akt. The homologous phosphorylation sites in Akt2 and Akt3 are S474 and S472, respectively.</p>

References:

1. Franke TF, et al., 1997, Cell 88: 435-7.
2. Burgering, BM and Coffey PJ, 1995, Nature 376: 599-602.
3. Franke TF, et al., 1995, Cell 81: 727-36.
4. Alessi DR, et al., 1996, EMBO J 15: 6541-51.
5. Sarbassov DD, et al. 2005, Science 307: 1098-101.
6. Jacinto E, et al., 2006, Cell 127: 125-37.
7. Cardone MH, et al., 1998, Science 282: 1318-21.
8. Brunet A, et al., 1999 Cell 96: 857-68.
9. Zimmermann S, and Moelling K, 1999 Science 286: 1741-4.
10. Cantley LC, and Neel BG, 1999 Proc Natl Acad Sci USA 96: 4240-5.
11. Vlahos CJ, et al., 1994, J Biol Chem 269: 5241-8.
12. Hajduch E, et al. 2001, FEBS Lett 492: 199-203.
13. Cross DA, et al., 1995 Nature 378: 785-9.
14. Diehl JA, et al., 1998 Genes Dev 12: 3499-511.
15. Gesbert F, et al., 2000, J Biol Chem 275: 39223-30.
16. Zhou BP, et al., 2001, Nat Cell Biol 3: 245-52.
17. Navé BT, et al., 1999 Biochem J 344: 427-31.
18. Inoki K, et al. 2002, Nat Cell Biol 4: 648-57.
19. Manning BD, et al., 2002, Mol Cell 10: 151-62.

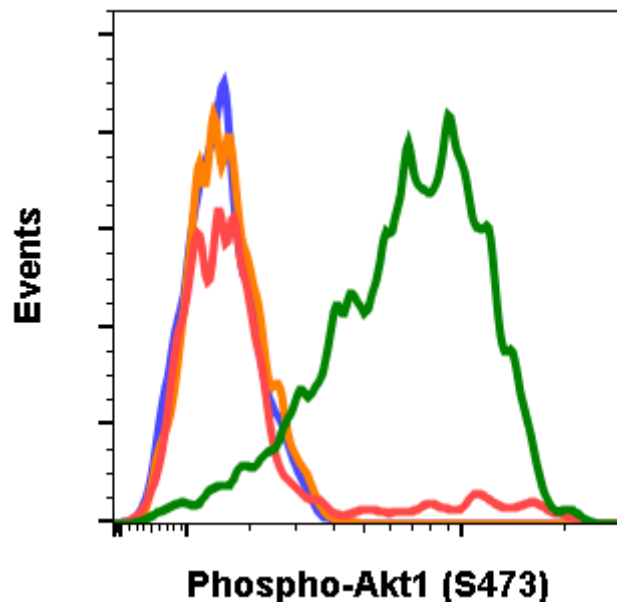

Flow cytometric analysis of C2C12 cells secondary antibody only negative control (blue) or 0.1 µg/mL of isotype control Cat. #2141 (orange) or untreated (red) or treated with staurosporine (green) using Phospho-Akt (Ser473) antibody AktS473-C7 at 0.1 µg/mL. Cat. #1131.

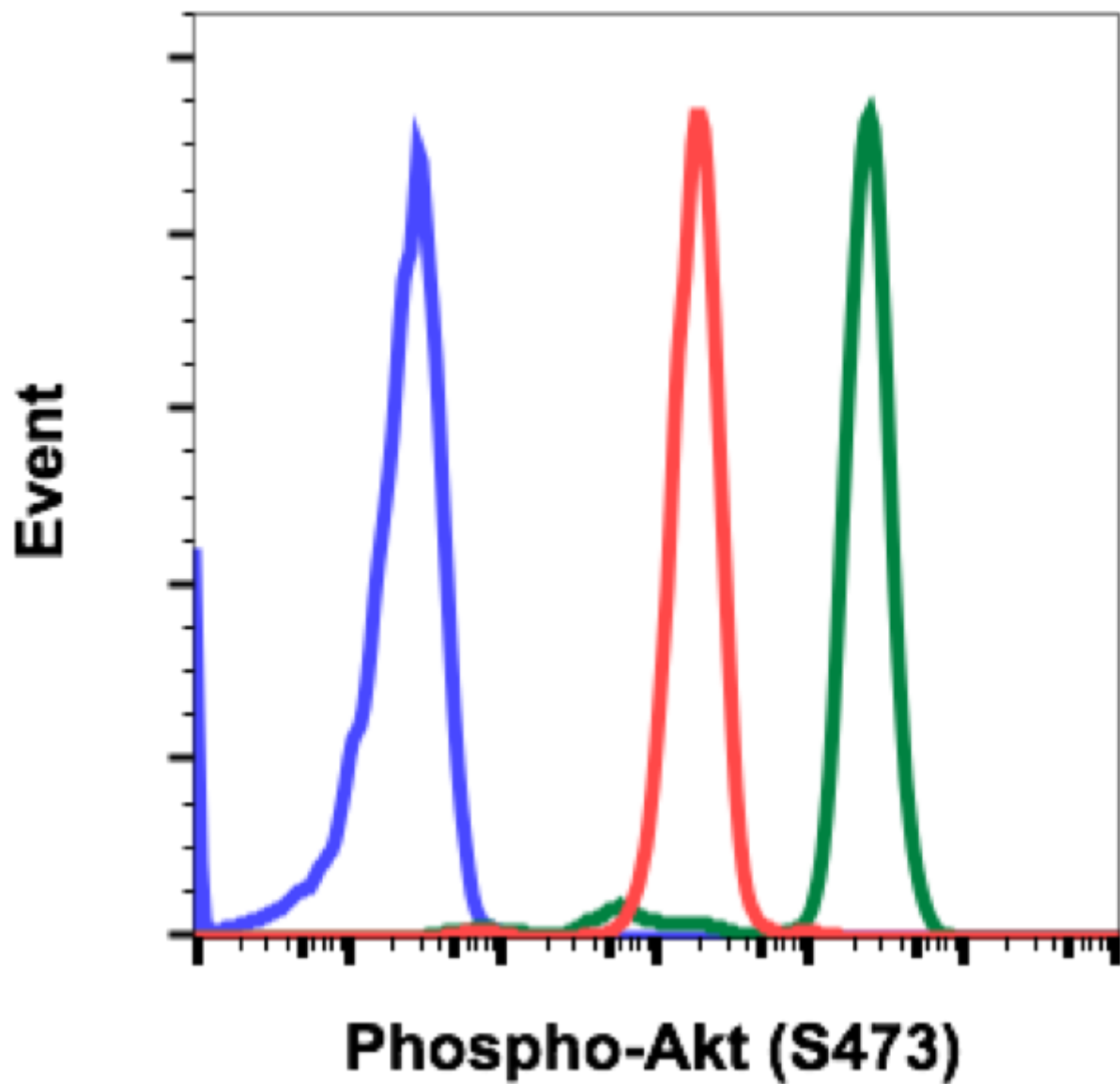

Flow cytometric analysis of U937 cells secondary antibody only negative control (blue) or treated with LY294002 (red) or with calyculin A (green) using Phospho-Akt (Ser473) antibody AktS473-C7. Cat. #1131.