

Phospho-Histone H3 (Ser10) (4B6) rabbit mAb FITC conjugate

www.abwizbio.com

Support: info@abwizbio.com

Order: sales@abwizbio.com

Catalog: #2063

Store at: 2-8°C

For Research Use Only. Not For Use In Diagnostic Procedures.

Applications

Flow Cytometry

Detection

N/A


Clonality

Monoclonal

Isotype

Rabbit IgGk

Format:	FITC
Cross Reactivity:	Predicted to work with mouse, rat and other homologues.
Formulation:	1X PBS, 0.09% NaN ₃ , 0.2% BSA
Preparation:	Protein A+G
Reactivity:	Human
Recommended Usage:	For flow cytometric staining, the suggested use of this reagent is 5 μ L per million cells or 5 μ L per 100 μ L of staining volume. It is recommended that the reagent be titrated for optimal performance for each application. See product image legends for additional information.
Immunogen:	A synthetic phospho-peptide corresponding to residues surrounding Ser10 of human phospho histone H3
Description:	Histones are highly conserved proteins that serve the core of nucleosomes, which serve to organize chromatin fiber for DNA packing. Histone H3 phosphorylation plays a major role in both transcriptional activation, which requires unpacking of the chromatin structure, and in chromosome packing during cell division. Histone H3 is phosphorylated at residues Ser10 and Ser28, and is acetylated at Lys14. Phosphorylation at Ser10 occurs during entry into mitosis prior to chromatin condensation, and phosphorylation at Ser28 follows a similar pattern. In response to EGF stimulation, it has been proposed that sequential Ser10 phosphorylation, then Lys14 acetylation occurs, causing a change in chromatin structure and gene activation.
References:	<ol style="list-style-type: none">1. Hans F and Dimitrov S. (2001) Oncogene. 20: 3021-3027.2. Cheung P, Tanner KG, Cheung WL, et al. (2000) Molecular Cell. 5: 905-915.

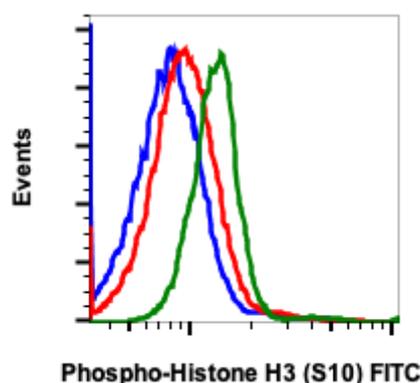

Flow cytometric analysis of HeLa cells unstained as negative control (blue) or untreated (red) or treated with nocodazole (green) and stained using Phospho-Histone H3 (Ser10) FITC conjugated antibody HisH3S10-4B6. Cat. #2063.