

Phospho-Stat5 (Tyr694) (G11) rabbit mAb SureLight 488 conjugate

www.abwizbio.com

Support: info@abwizbio.com

Order: sales@abwizbio.com

#2330

Store at: 2-8°C

For Research Use Only. Not For Use In Diagnostic Procedures.

Applications	Detection	Clonality	Isotype
Flow Cytometry	N/A	Monoclonal	Rabbit IgGk

Format: SureLight 488

Cross Reactivity: Predicted to work with mouse, rat and other homologues.

Formulation: 1X PBS, 0.09% NaN₃, 0.2% BSA

Preparation: Protein A+G

Reactivity: Human, Mouse

Recommended


Usage: For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells or 5 µL per 100 µL of staining volume. It is recommended that the reagent be titrated for optimal performance for each application.

Immunogen: A synthetic phospho-peptide corresponding to residues surrounding Tyr695 of human phospho Stat5

Description: Stat5 activation occurs in response to many ligands including prolactin, IL-2, growth hormone, and GM-CSF. Tyr694 phosphorylation is obligatory activation of Stat5 (1,2), and is mediated by Src upon erythropoietin stimulation (3). Phospho Stat5 is constitutively active in some leukemic cell types (4), and phospho Stat5 is found in some endothelial cells when treated with IL-3, suggesting its involvement in cell motility and angiogenesis (5). Stat5 has been shown to be encoded by two separate genes, Stat5a and Stat5b, which share over 90% amino acid sequence identity. In different cell types, Stat5a and Stat5b are independently regulated and activated. For example, interferon treatment predominantly activates Stat5a in U937 cells and Stat5b in HeLa cells (6).

References:

1. Gouilleux, F. et al. (1994) EMBO J. 13:4361-4369.
2. Wakao, H. et al. (1994) EMBO J. 13:2182-2191.
3. Okutani, Y. et al. (2001) Oncogene. 20:6643-6650.
4. Demoulin, J.B. et al. (1999) J. Biol. Chem. 274:25855-25861.
5. Dentelli, P. et al. (1999) J. Immunol. 163:2151-2159.
6. Meinke, A. et al. (1996) Mol. Cell. Biol. 16:6937-6944.

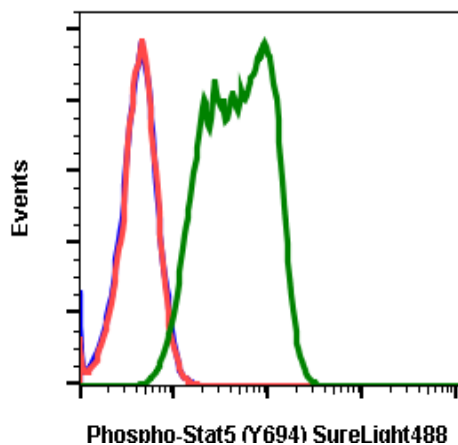

Flow cytometric analysis of NIH3T3 cells unstained and treated with imatinib as negative control (blue) or stained and treated with imatinib (red) or treated with pervanadate (green) using Phospho-Stat5 (Tyr694) antibody Stat5Y694-G11 SureLight488 conjugate. Cat. #2330.