

Phospho-Btk (Tyr551) (G12) rabbit mAb SureLight488 conjugate

www.abwizbio.com

Support: info@abwizbio.com

Order: sales@abwizbio.com

Catalog: #2345

Store at: 2-8°C

For Research Use Only. Not For Use In Diagnostic Procedures.

Applications

Flow Cytometry

Detection

N/A


Clonality

Monoclonal

Isotype

Rabbit IgGk

Format:	SureLight 488
Cross Reactivity:	Predicted to work with mouse, rat and other homologues.
Formulation:	1X PBS, 0.09% NaN ₃
Preparation:	Protein A+G
Reactivity:	Human
Recommended Usage:	For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells or 5 µL per 100 µL of staining volume. It is recommended that the reagent be titrated for optimal performance for each application.
Immunogen:	A synthetic phosphor-peptide corresponding to residues surrounding Tyr551 of human phospho Btk
Description:	Btk is a major node in the B-cell receptor signaling pathway, where it regulates B cell maturation, activation, survival, differentiation, and proliferation. Btk is activated by Src family kinases, including Lyn, which phosphorylates Btk at Tyr551. Upon phosphorylation at this site, Btk is recruited to the plasma membrane where autophosphorylation at Tyr223 occurs. The Btk signaling pathway is a major target of small molecule inhibitors for B-cell lymphoma, autoimmune diseases, and non-Hodgkin's lymphomas. These inhibitors either form a covalent bond at Cys481 in the ATP-binding site or serve as reversible inhibitors that bind the SH3 pocket and stabilize inactive Btk.
References:	Liang, C, Tian D, Ren X, et al. (2018) European Journal of Medicinal Chemistry. 151: 315-326.

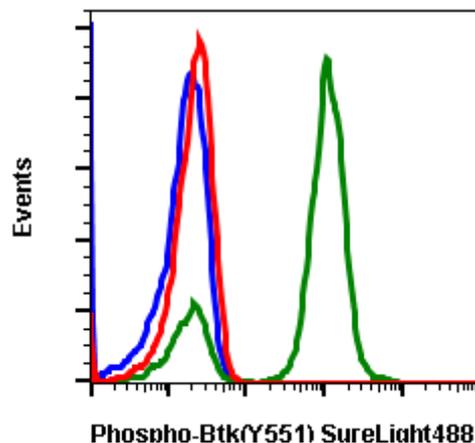

Flow cytometric analysis of Ramos cells unstained and untreated as negative control (blue) or stained untreated (red) or treated with pervanadate (green) using Phospho-Btk (Tyr551) antibody BtkY551-G12 SureLight488 conjugate. Cat. #2345.