

Phospho-Shp1 (Tyr536) (2A7) rabbit mAb FITC conjugate

www.abwizbio.com

Support: info@abwizbio.com

Order: sales@abwizbio.com

Catalog: #2353

Store at: 2-8°C

For Research Use Only. Not For Use In Diagnostic Procedures.

Applications

Flow Cytometry

Detection

N/A


Clonality

Monoclonal

Isotype

Rabbit IgGk

Format:	FITC
Cross Reactivity:	Predicted to work with mouse, rat and other homologues.
Formulation:	1X PBS, 0.09% NaN ₃ , 0.2% BSA
Preparation:	Protein A+G
Reactivity:	Human
Recommended Usage:	For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells or 5 µL per 100 µL of staining volume. It is recommended that the reagent be titrated for optimal performance for each application. See product image legends for additional information.
Immunogen:	A synthetic phospho-peptide corresponding to residues surrounding Tyr536 of human phospho Shp1
Description:	Src homology-2 domain containing phosphatase-1 (Shp1) typically functions as a negative regulator for signal transduction, interacting with Zap70, CD5, CD3ε, and IL-2R in T cells. Shp1 suppresses cellular signaling by dephosphorylating target proteins. Shp1 is expressed predominantly in hematopoietic cells, where it can play a role in malignant transformation into cancer cells. Methylation of the Shp1 promoter has been shown to silence Shp1 expression and promote the development of cancer. Shp1 negatively regulates PDGF and insulin tyrosine receptor pathways, and negatively regulates Lck in the TCR signaling pathway to limit TCR responsiveness. Shp1 activity is prevented by TAOK3 in T cells to allow for T cell activation and differentiation.
References:	<ol style="list-style-type: none">1. Oka T, Ouchida M, Koyama M, et al. (2002) Cancer Research. 62:6390-6394.2. Gerald P, Hiraoka-Yamamoto J, Matsumoto M, et al., (2009) Nature Medicine. 15:1298-1306.3. Ormonde JVS, Li Z, Stegen C, and Madrenas J. (2018) Journal of Immunology 201:3431-3442.

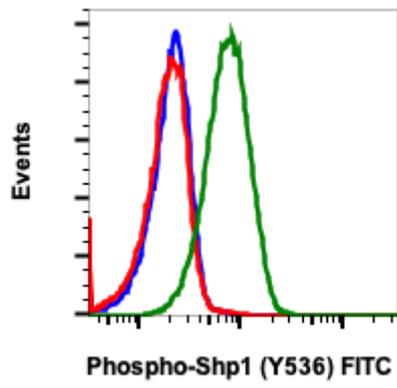

Flow cytometric analysis of Ramos cells unstained as negative control (blue) or untreated (red) or treated with pervanadate (green) and stained using Phospho-Shp1 (Tyr536) antibody Shp1Y536-2A7 FITC conjugate. Cat#2353.