

# Phospho-SAPK/JNK (Thr183/Tyr185) (A11) rabbit mAb PE conjugate

[www.abwizbio.com](http://www.abwizbio.com)

Support: [info@abwizbio.com](mailto:info@abwizbio.com)

Order: [sales@abwizbio.com](mailto:sales@abwizbio.com)


Catalog: #2357

Store at: 2-8°C

For Research Use Only. Not For Use In Diagnostic Procedures.

Applications	Detection	Clonality	Isotype
Flow Cytometry	N/A	Monoclonal	Rabbit IgGk

<b>Format:</b>	PE
<b>Cross Reactivity:</b>	Predicted to work with mouse, rat and other homologues.
<b>Formulation:</b>	1X PBS, 0.09% NaN <sub>3</sub> , 0.2% BSA
<b>Preparation:</b>	Protein A+G
<b>Reactivity:</b>	Human
<b>Recommended Usage:</b>	For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells or 5 µL per 100 µL of staining volume. It is recommended that the reagent be titrated for optimal performance for each application. See product image legends for additional information.
<b>Immunogen:</b>	A synthetic phospho-peptide corresponding to residues surrounding Thr183/Tyr185 of human phospho SAPK/JNK
<b>Description:</b>	The SAPK/JNK pathway initiates apoptosis upon exposure to radiation, UV exposure, heat shock, oxidative stress, and other stressors. Upon exposure to environmental stress, the SAPK/JNK signaling pathway sequentially activates the proteins MEKK1, SEK1, SAPK, and c-Jun. Upstream activators of the SAPK/JNK cascade include ceramide, small GTP-binding proteins such as Rac1 and Cdc42Hs, Ask1, and caspases. MKK7 is also a major and direct SAPK/JNK activator in the TNF $\alpha$ or environmental stress signaling pathways, where its kinase activity directly phosphorylates SAPK/JNK. This relationship between MKK7 and SAPK appears to be evolutionarily conserved, as it is preserved in their Drosophila homologues, Hep and DJNK, respectively.
<b>References:</b>	Verheij M, Bose R, Lin XH, et al. (1996) Nature. 380: 75-79. Verheij M, Ruiter GA, Zerp SF, et al. (1998) Radiotherapy and Oncology. 47: 225-232. Moriguchi T, Toyoshima F, Masuyama N, et al. (1997) The EMBO J 16: 7045-7053.

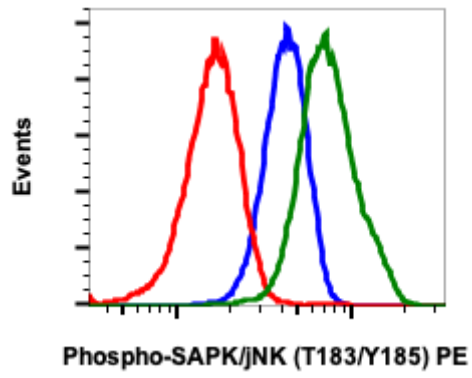

Flow cytometric analysis of 293T cells untreated and unstained as negative control (blue) or untreated (red) or with UV+TPA (green) and stained using Phospho-SAPK/JNK (Thr183/Tyr185) antibody SAPKT183Y185-A11 PE conjugate. Cat. #2357.