

Phospho-NFκB p65 (Ser529) (H3) rabbit mAb PE conjugate

www.abwizbio.com

Support: info@abwizbio.com

Order: sales@abwizbio.com

Catalog: #2407

Store at: 2-8°C

For Research Use Only. Not For Use In Diagnostic Procedures.

Applications

Flow Cytometry

Detection

N/A


Clonality

Monoclonal

Isotype

Rabbit IgGκ

Format:	PE
Cross Reactivity:	Predicted to work with mouse, rat and other homologues.
Formulation:	1X PBS, 0.09% NaN ₃ , 0.2% BSA
Preparation:	Protein A+G
Reactivity:	Human, Mouse
Recommended Usage:	For flow cytometric staining, the suggested use of this reagent is 5 μL per million cells or 5 μL per 100 μL of staining volume. It is recommended that the reagent be titrated for optimal performance for each application. See product image legends for additional information.
Immunogen:	A synthetic phospho-peptide corresponding to residues surrounding Ser529 of human phospho-NFKB p65
Description:	The nuclear factor κB (NFκB)/Rel family of transcription factors play a pivotal role in inflammatory and immune responses (1,2). NF-kappa-B is present in almost all cell types and is involved in many biological processes including immunity, inflammation, cell growth and differentiation, apoptosis, and tumorigenesis. NFκB is a homo- or heterodimeric complex formed by the Rel-like domain-containing proteins RELA/p65, RELB, NFκB1/p105, NFκB1/p50, REL and NFκB2/p52. The dimers bind at κB sites in the target gene DNA. Individual dimers have distinct preferences for different κB sites and can act as either transcriptional activators or repressors. NFκB Ser536 phosphorylation stimulates Lys310 acetylation and interaction of phospho NFκB with CBP. Acetylated/phospho NFκB induces enhanced transcriptional activity.
References:	<ol style="list-style-type: none">1. Baeuerle, P.A. and Henkel, T. (1994) Annu. Rev. Immunol. 12:141-179.2. Baeuerle, P.A. and Baltimore, D. (1996) Cell. 87:13-20.

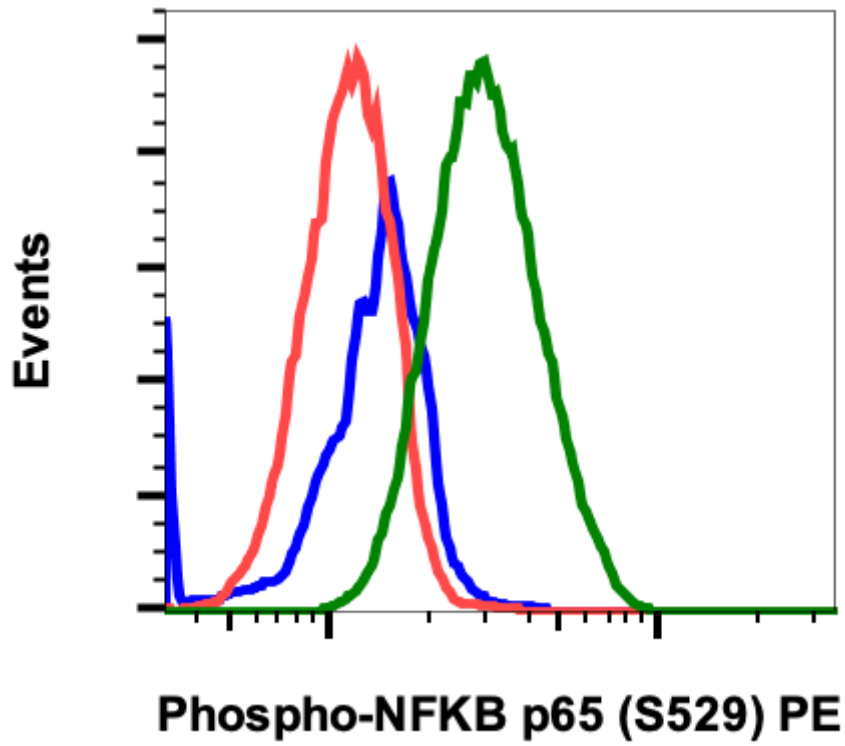

Flow cytometric analysis of HeLa cells untreated and unstained as negative control (blue) or untreated (red) or treated with TNF α + calyculin A (green) and stained using Phospho-NFkB p65 (Ser529) antibodyPE conjugate, NFkBp65S529-H3. Cat. #2407.