

Phospho-MARCKS (Ser167/170) (C9) rabbit mAb FITC conjugate

www.abwizbio.com

Support: info@abwizbio.com

Order: sales@abwizbio.com


Catalog: #2448

Store at: 2-8°C

For Research Use Only. Not For Use In Diagnostic Procedures.

Applications	Detection	Clonality	Isotype
Flow Cytometry	N/A	Monoclonal	Rabbit IgGκ

Format:	FITC
Cross Reactivity:	Predicted to work with mouse, rat and other homologues.
Formulation:	1X PBS, 0.09% NaN ₃ , 0.2% BSA
Preparation:	Protein A+G
Reactivity:	Human,Mouse,Rat
Recommended Usage:	For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells or 5 µL per 100 µL of staining volume. It is recommended that the reagent be titrated for optimal performance for each application. See product image legends for additional information.
Immunogen:	A synthetic phospho-peptide corresponding to residues surrounding Ser167/170 of human phospho MARCKS
Description:	MARCKS (myristoylated alanine-rich C kinase substrate) is a major PKC substrate expressed in all eukaryotic cells(1,2). It binds to and cross-links actin filaments to serve as a bridge between Ca ²⁺ /calmodulin and PKC signaling and attenuates phosphatidylinositol 4,5-bisphosphate plasma membrane signaling (3). MARCKS is involved with cell mobility, phagocytosis, membrane traffic, cell adhesion, and mitogenesis. Ser159, 163, 167 and 170 of MARCKS are phosphorylated by PKC in response to cell growth and cellular stress (4). MARCKs phosphorylation is believed to induce its translocation from plasma membrane to cytoplasm.
References:	<ol style="list-style-type: none">1. El Amri M et al., (2018) J Biomed Sci 25(1):43. doi: 10.1186/s12929-018-0445-2. Aderem A. (1992) Cell 71:713-6.3. Hartwig JH, et al., (1992) Nature 356:618-22.4. Bhat NR. et al., (1991) J Neurosci Res 30: 447-54

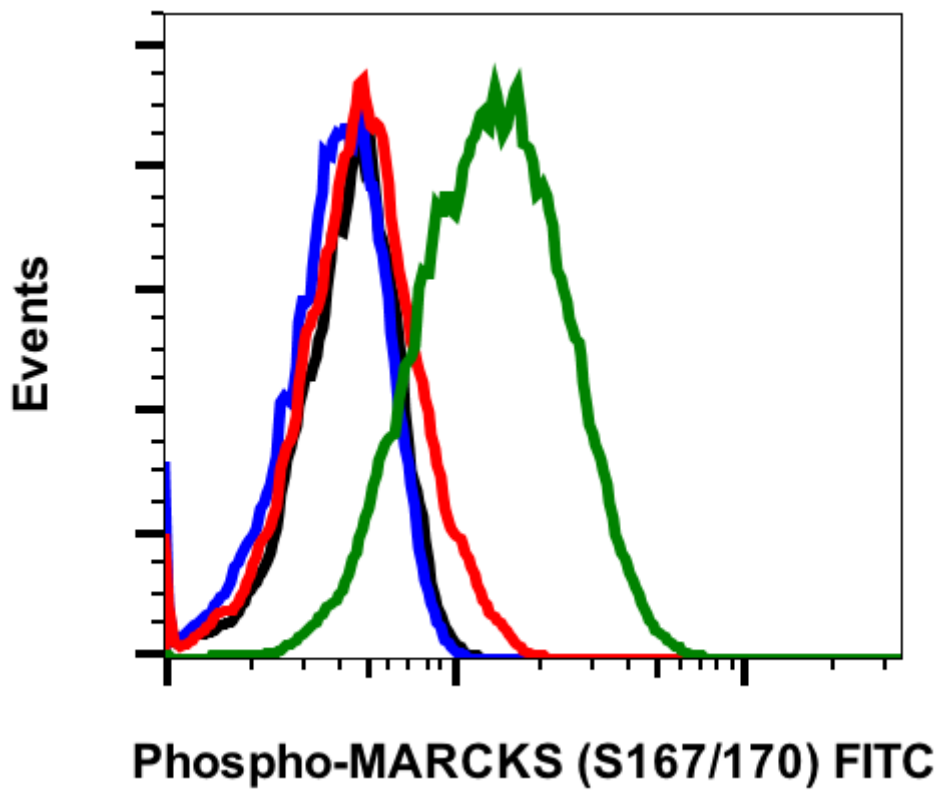

Flow cytometric analysis of C6 cells treated with staurosporine (red) or treated with UV+TPA (green) using Phospho-MARCKS (Ser167/170) (C9) Rabbit mAb (FITC Conjugate) MARKSS167170-C9 #2448, or concentration-matched Rabbit (G9) mAb IgG Isotype Control (FITC Conjugate) #2143 for cells treated with staurosporine (black) or treated with UV TPA (blue).