

SARS-CoV-2 Alpha B.1.1.7 Variant Recombinant Spike Trimer His Tag

www.abwizbio.com

Support: info@abwizbio.com

Order: sales@abwizbio.com

#2581

Store at: 2-8°C

For Research Use Only. Not For Use In Diagnostic Procedures.

Applications	Detection	Clonality	Isotype
ELISA	Anti-SARS-CoV-2 mAb	SARS-CoV-2	SARS-CoV-2

Format: His tag

Cross Reactivity: Predicted to work with mouse, rat and other homologues.

Formulation: 1X PBS, 0.025% NaN₃

Preparation: His tag purification

Reactivity: Other

Recommended


Usage: Control antigen for SARS-CoV-2 (COVID-19) diagnostic assays. Contains C-terminal His tag used for purification.

Immunogen: N/A

Description: N/A


References:

1. Walls A.C. et al., 2019 Cell, 176: 1026-1039.
2. Tang T. et al., 2020, Antiviral Res., 178:10479.
3. Jiang S et al., 2022, Trends Immunol, 41:355-359.
4. Bosch BJ, et al., 2003, J Virol 77:8801-8811.

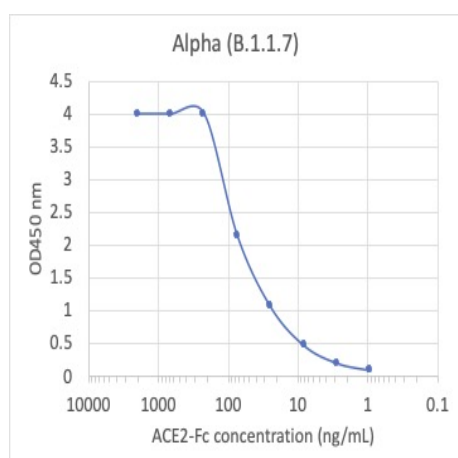

Concentration?response curves for binding of CoV2 spike protein to human ACE2 in cell-free ELISA-type assays. Microtiter wells were coated with 100 μ L of each spike trimer at 2 μ g/mL in PBS at 4°C overnight. The wells were washed with PBS and blocked with 200 μ L of 1% BSA/PBS. ACE2-Fc was serially diluted from 2 μ g/mL in 1% BSA/PBS. The blocker was discarded, and the wells were incubated with 100 μ L of serially diluted ACE2-Fc at 37°C for 1 hour. The wells were washed with PBS and the bound ACE2-Fc was detected with 100 μ L of Peroxidase AffiniPure Goat Anti-Human IgG, Fc γ fragment specific (Jackson Immuno Research 109-035-098)(1:5,000 in 1% BSA/PBS) at 37°C for 1 hour. The wells were washed with PBS and the wells were developed with 100 μ L of MB/E Ultra Sensitive, Blue, Horseradish Peroxidase Substrate (EMD Millipore ES022-500ml) at RT for 5 min. The reaction was stopped with 100 μ L of 0.6N H₂SO₄ and the signals were read at 450 nm using a plate reader (Biotek).

**SARS-CoV2 Alpha (B.1.1.7)
Recombinant Spike
Trimer Protein His Tag**

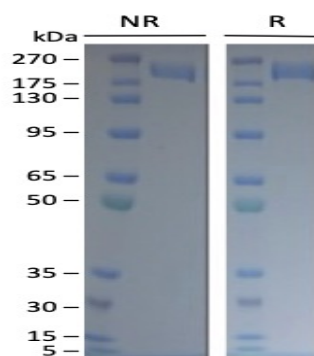

SARS-CoV-2 recombinant spike trimer His Tag protein (Alpha| B.1.1.7) on SDS-PAGE under non-reducing (NR) and reducing (R) conditions. The gel was stained with Coomassie G-250 per manufacturer instructions. The purity of the protein is greater than 95% and it is between 175 and 270 kDa.