

# SARS-CoV-2 BA.2.75 Omicron Variant Recombinant Spike Trimer His Tag

[www.abwizbio.com](http://www.abwizbio.com)  
Support: [info@abwizbio.com](mailto:info@abwizbio.com)  
Order: [sales@abwizbio.com](mailto:sales@abwizbio.com)

Catalog: #2664


Store at: 2-8°C

*For Research Use Only. Not For Use In Diagnostic Procedures.*

Applications	Detection	Clonality	Isotype
ELISA	Anti-SARS-CoV-2 mAb	SARS-CoV-2	SARS-CoV-2


Format:	His tag
Cross Reactivity:	Predicted to work with mouse, rat and other homologues.
Formulation:	1X PBS, 0.025% NaN3
Preparation:	His tag purification
Reactivity:	Other
Recommended Usage:	Control antigen for SARS-CoV-2 (COVID-19) diagnostic assays. Contains C-terminal His tag used for purification.
Immunogen:	N/A
Description:	
References:	1. Walls A.C. et al., 2019 Cell, 176: 1026-1039. 2. Tang T. et al., 2020, Antiviral Res., 178:10479. 3. Jiang S et al., 2022, Trends Immunol, 41:355-359. 4. Bosch BJ, et al., 2003, J Virol 77:8801-8811.

### BA.2.75 spike trimer binding to ACE2

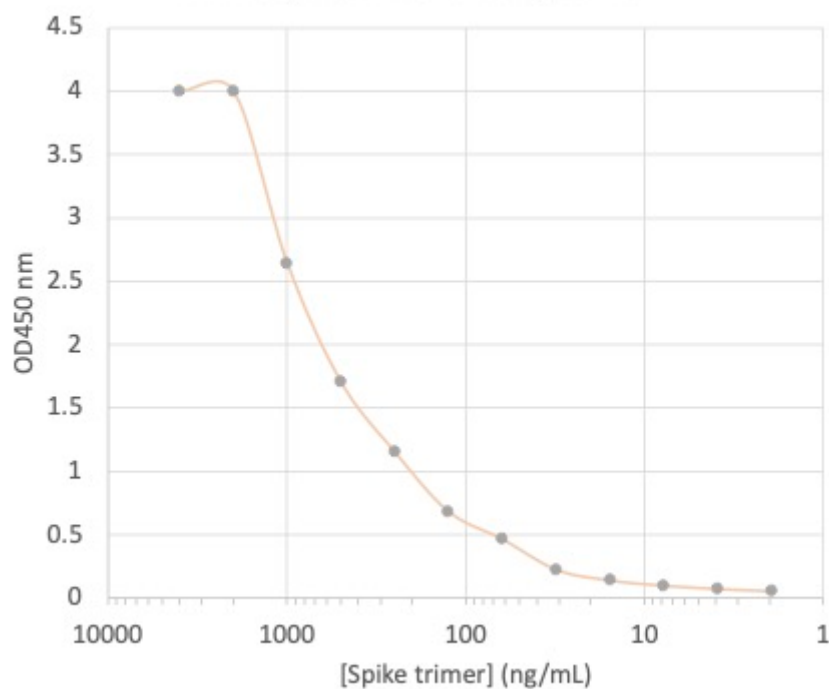

Concentration-response curves for binding of CoV2 spike protein to human ACE2 in cell-free ELISA-type assays. Microtiter wells were coated with 100  $\mu$ L of ACE2-Fc at 2  $\mu$ g/mL in PBS at 4°C overnight. The wells were washed with PBS and blocked with 200  $\mu$ L of 1% BSA/PBS at RT for 1 hour. Spike trimer was serially diluted from 10  $\mu$ g/mL in 1% BSA/PBS. The blocker was discarded, and the wells were incubated with 100  $\mu$ L of serially diluted spike trimer at RT for 1 hour. The wells were washed with PBS and the bound spike trimer was detected with 100  $\mu$ L of Peroxidase AffiniPure Rabbit Anti-His Tag (Jackson Immuno Research 300-035-240)(1:5,000 in 1% BSA/PBS) at RT for 1 hour. The wells were washed with PBS and the wells were developed with 100  $\mu$ L of MB/E Ultra-Sensitive, Blue, Horseradish Peroxidase Substrate (EMD Millipore ES022-500ml) at RT for 5 min. The reaction was stopped with 100  $\mu$ L of 0.6N H<sub>2</sub>SO<sub>4</sub> and the signals were read at 450 nm using a plate reader (Biotek).

### SARS-CoV-2 BA.2.75 Spike Timer Protein

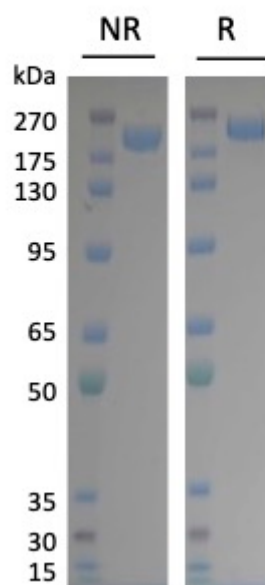

SARS-CoV-2 Spike Trimer, His Tag (BA.2.75/Omicron) on SDS-PAGE under non-reducing (NR) and reducing (R) conditions. The gel was stained with Coomassie G-250 per manufacturer instructions. The purity of the protein is greater than 95% and it is between 175 and 270 kDa.