

SARS-CoV-2 BA.2.75.2 Omicron Variant Recombinant Spike Trimer His Tag

www.abwizbio.com
Support: info@abwizbio.com
Order: sales@abwizbio.com

Catalog: #2676


Store at: 2-8°C

For Research Use Only. Not For Use In Diagnostic Procedures.

Applications	Detection	Clonality	Isotype
ELISA	Anti-SARS-CoV-2 mAb	SARS-CoV-2	SARS-CoV-2

Format:	His tag
Cross Reactivity:	Predicted to work with mouse, rat and other homologues.
Formulation:	1X PBS, 0.025% NaN ₃
Preparation:	His tag purification
Reactivity:	Other
Recommended Usage:	Control antigen for SARS-CoV-2 (COVID-19) diagnostic assays. Contains C-terminal His tag used for purification.
Immunogen:	N/A
Description:	
References:	<ol style="list-style-type: none">1. Walls A.C. et al., 2019 Cell, 176: 1026-1039.2. Tang T. et al., 2020, Antiviral Res., 178:10479.3. Jiang S et al., 2022, Trends Immunol, 41:355-359.4. Bosch BJ, et al., 2003, J Virol 77:8801-8811.


BA.2.75.2 spike trimer binding to ACE2

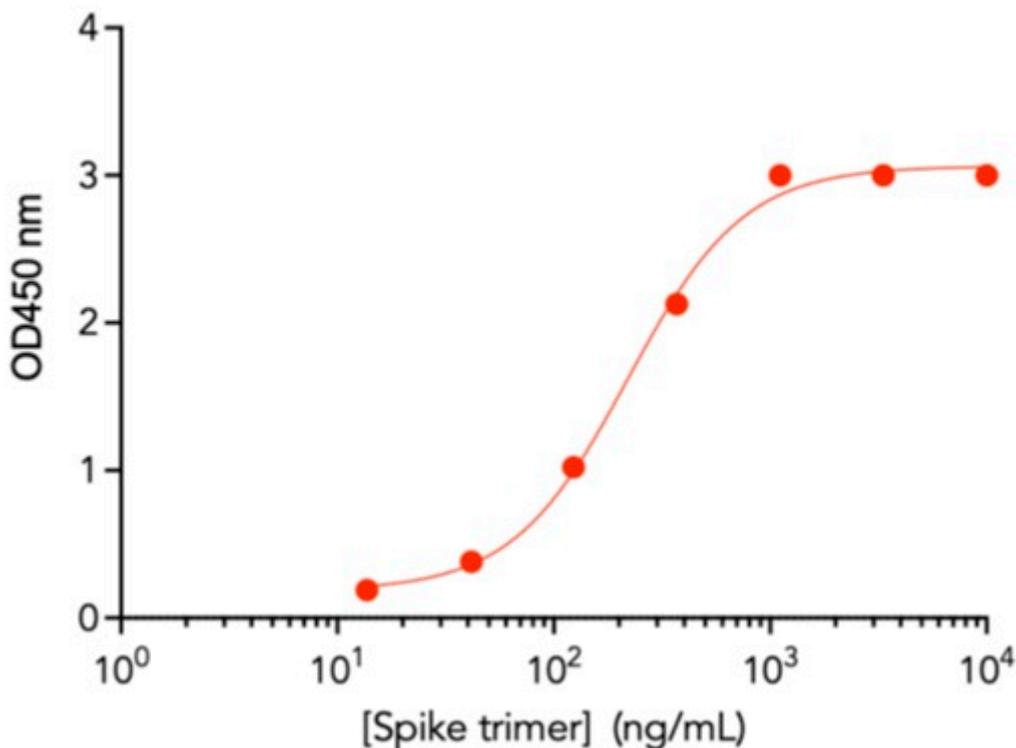

Concentration-response curves for binding of CoV2 spike protein to human ACE2 in cell-free ELISA-type assays. Microtiter wells were coated with 100 μ L of ACE2-Fc at 2 μ g/mL in PBS at 4°C overnight. The wells were washed with PBS and blocked with 200 μ L of 1% BSA/PBS at RT for 1 hour. Spike trimer was serially diluted from 10 μ g/mL in 1% BSA/PBS. The blocker was discarded, and the wells were incubated with 100 μ L of serially diluted spike trimer at RT for 1 hour. The wells were washed with PBS and the bound spike trimer was detected with 100 μ L of Peroxidase AffiniPure Rabbit Anti-His Tag (Jackson Immuno Research 300-035-240)(1:5,000 in 1% BSA/PBS) at RT for 1 hour. The wells were washed with PBS and the wells were developed with 100 μ L of MB/E Ultra-Sensitive, Blue, Horseradish Peroxidase Substrate (EMD Millipore ES022-500ml) at RT for 5 min. The reaction was stopped with 100 μ L of 0.6N H₂SO₄ and the signals were read at 450 nm using a plate reader (Biotek).

SARS-CoV-2 Omicron BA.2.75.2 Variant

Recombinant Spike Trimer His Tag

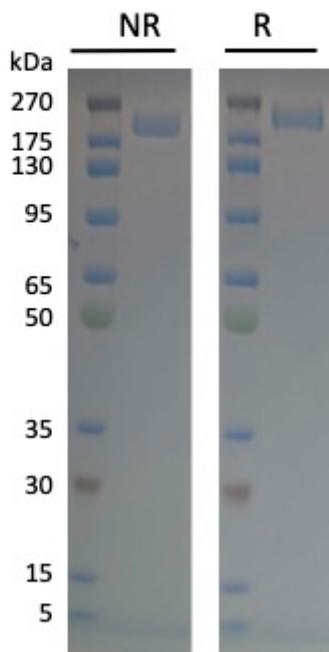

SARS-CoV-2 Spike Trimer, His Tag (BA.2.75.2) on SDS-PAGE under non-reducing (NR) and reducing (R) conditions. The gel was stained with Coomassie G-250 per manufacturer instructions. The purity of the protein is greater than 95% and it is between 175 and 270 kDa.